

**GORE® TAG®**


Conformable Thoracic Stent Graft  
with ACTIVE CONTROL System



# Expanded sizes, proven performance

Offer more patients the treatment you trust  
with our large-diameter tapers.

Optimize oversizing for a broad range of tapered aortas

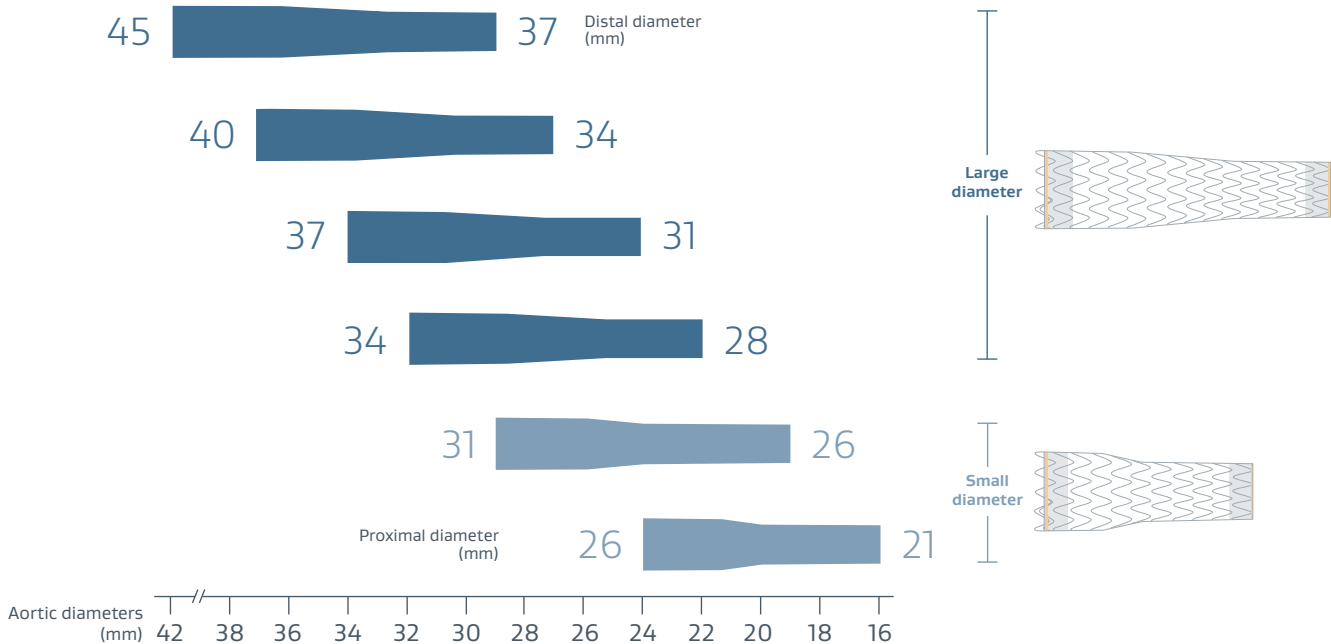

*Together, improving life*



# A design against dSINE

Tapered designs may reduce distal stent graft-induced new entry (dSINE) tears associated with oversizing in a narrow true lumen.<sup>1-3</sup>

dSINE is linked to increased reinterventions and rupture risk<sup>1-3</sup>

## Tapered devices


Catalogue number*	Labeled diameter (mm)	Intended aortic diameter (mm)	Device length (cm)	Device profile (Fr)	Oversizing range (%)	Partially uncovered stent length (mm)
TGM262110	26 x 21	19.5-24/16-19.5	10	20	8-33	4
TGMR312610	31 x 26	24-29/19.5-24	10	20	7-33	4
TGM342815	34 x 28	27-32/22-26	15	22	6-27	5
TGMR373115	37 x 31	29-34/24-29	15	22	7-29	5
TGMR403415	40 x 34	31-37/27-32	15	22	6-29	6
TGM453715	45 x 37	34-42/29-34	15	24	7-32	6.5

\* Some sizes are not available in all regions.

## References

1. Canaud L, Gandet T, Sfeir J, Ozdemir BA, Solovei L, Alric P. Risk factors for distal stent graft-induced new entry tear after endovascular repair of thoracic aortic dissection. *Journal of Vascular Surgery* 2019;69(5):1610-1614. OPEN ACCESS: <https://www.sciencedirect.com/science/article/pii/S0741521418322560>
2. János RA, Tsagakis K, Bettin M, et al. Thoracic aortic aneurysm expansion due to late distal stent graft-induced new entry. *Catheterization & Cardiovascular Interventions* 2015;85(2):E43-E53.
3. Lortz J, Leinburger F, Tsagakis K, et al. Distal stent graft induced new entry: risk factors in acute and chronic type B aortic dissections. *European Journal of Vascular & Endovascular Surgery* 2019;58(6):822-830. OPEN ACCESS <https://www.sciencedirect.com/science/article/pii/S1078588419302928>

 Consult Instructions for Use  
[eifu.goremedical.com](http://eifu.goremedical.com)

**INDICATIONS FOR USE IN THE U.S.:** The GORE® TAG® Conformable Thoracic Stent Graft is intended for endovascular repair of all lesions of the descending thoracic aorta, including: isolated lesions in patients who have appropriate anatomy, including: adequate iliac/femoral access, aortic inner diameter in the range of 16-42 mm, ≥ 20 mm non-aneurysmal aorta proximal and distal to the lesion; Type B dissections in patients who have appropriate anatomy, including: adequate iliac/femoral access, ≥ 20 mm landing zone proximal to the primary entry tear; proximal extent of the landing zone must not be dissected, diameter at proximal extent of proximal landing zone in the range of 16-42 mm. **CONTRAINDICATIONS:** Patients with known sensitivities or allergies to the device materials; patients who have a condition that threatens to infect the graft. Refer to *Instructions for Use* at [eifu.goremedical.com](http://eifu.goremedical.com) for a complete description of all applicable indications, warnings, precautions and contraindications for the markets where this product is available. 

Products listed may not be available in all markets.

GORE, *Together, improving life*, ACTIVE CONTROL, TAG and designs are trademarks of W. L. Gore & Associates.

© 2025 W. L. Gore & Associates, Inc. 24AR1217-EN03 FEBRUARY 2025

W. L. Gore & Associates, Inc.  
[goremedical.com](http://goremedical.com)

Asia Pacific +65 67332882 Australia/New Zealand 1800 680 424 Europe 00800 6334 4673  
United States Flagstaff, AZ 86004 800 437 8181 928 779 2771

