

PRECISE HANDLING

Without sacrificing tissue tolerance



GORE-TEX® Suture is a microporous, nonabsorbable monofilament made of expanded polytetrafluoroethylene (ePTFE).

■ **Tissue Tolerance**

The biocompatibility and inertness of ePTFE allows for the GORE-TEX® Suture to remain in the oral environment for as long as two to four weeks.

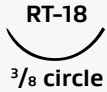
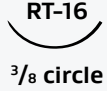

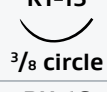
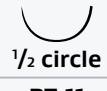
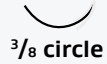
■ **Handling**

Low friction and the smooth, supple nature of the GORE-TEX® Suture allow for precise handling and provide flexibility in the positioning of a square knot.

■ **Nonwicking**

The monofilament GORE-TEX® Suture is not subject to bacterial wicking sometimes associated with multifilament sutures.

Configurations

Catalog Number	Needle Size	Thread Size	Description
P4K13E	 RT-18 3/8 circle	CV-4	A CV-4 suture with an 18 mm reverse-cutting needle. This can be used as an alternative to the P5K23 where a stronger suture is desired.
P5K17E	 RT-16 3/8 circle	CV-5	A CV-5 Suture with a 16mm reverse-cutting needle. Many clinicians select this suture for dental implant procedures or flap procedures in which they would prefer to leave the sutures in place for extended periods of time.
P5K23E	 RT-18 3/8 circle	CV-5	A CV-5 suture with an 18 mm reverse-cutting needle. This can be used as an alternative to the P5K17 where a longer needle is preferred.
P6K23E	 RT-13 3/8 circle	CV-6	A CV-6 suture with a 13 mm reverse-cutting needle. This is a finer suture with a smaller needle for delicate procedures such as gingival grafts or mucosal suturing.
P6K25E	 RH-16 1/2 circle	CV-6	A CV-6 suture with a 16 mm reverse-cutting needle. This can be used as an alternative to the P6K23 where a 1/2 circle needle is preferred.
P7K13E	 RT-11 3/8 circle	CV-7	A CV-7 suture with an 11 mm reverse-cutting needle. This is a finer suture with a smaller needle for delicate procedures such as gingival grafts or mucosal suturing.

Actual Size

* Not USP/Not E.P.

 Consult Instructions
for Use
eifu.goremedical.com

Refer to *Instructions for Use* at eifu.goremedical.com for a complete description of all applicable indications, warnings, precautions and contraindications for the markets where this product is available. ®/©

Products listed may not be available in all markets.

GORE, GORE-TEX, *Together, improving life* and designs are trademarks of W. L. Gore & Associates.
© 1996, 1998, 2000, 2001, 2005, 2007, 2008, 2012, 2014, 2023 W. L. Gore & Associates GmbH 231227957-EN OCTOBER 2023

W. L. Gore & Associates, Inc.
goremedical.com

Asia Pacific +65 6733 2882 Australia/New Zealand 1800 680 424 Europe 00800 6334 4673
United States Flagstaff, AZ 86004 800 437 8181 928 779 2771

