

# TIME TO RETHINK TIPS



GORE® VIATORR® TIPS Endoprosthesis with Controlled Expansion combines the legacy of proven patency<sup>1</sup> with diameter control

- **Control the diameter** — Designed to reach a targeted portal pressure gradient
- **Lasting diameter control\*** — Size and set the diameter to stay
- **Engineered for flexibility** — Conformability to tortuous anatomy
- **Single device** — Size to a diameter of 8 / 9 / 10 mm during implantation

\* Based on benchtop data on file. Less than 0.25 mm increase in diameter (diameter expansion) demonstrated by a simulated 10 year period at physiologic portal pressures.

1. Based on GORE® VIATORR® TIPS Endoprosthesis: Bureau C, Pagan JCG, Layrargues GP, et al. Patency of stents covered with polytetrafluoroethylene in patients treated by transjugular intrahepatic portosystemic shunts: long term results of a randomized multicentre study. *Liver International* 2007;27(6):742–747.

*Together, improving life*








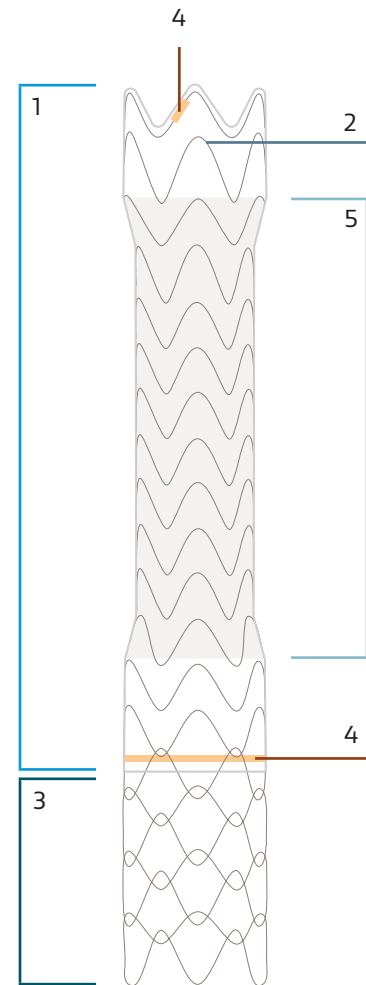
## Device specifications

Catalogue number	Endoprosthesis internal diameter (mm)*	Graft-lined length (cm)	Unlined length (cm)	Maximum guidewire diameter (in)	Hemostatic introducer sheath (Fr)	Maximum dilatation balloon diameter (mm)*
PTB8104275	8–10	4	2	≤ 0.035	10	10
PTB8105275	8–10	5	2	≤ 0.035	10	10
PTB8106275	8–10	6	2	≤ 0.035	10	10
PTB8107275	8–10	7	2	≤ 0.035	10	10
PTB8108275	8–10	8	2	≤ 0.035	10	10

\* The selected balloon diameter may be 8, 9, or 10 mm, and should not exceed 10 mm. A balloon which can reach an inflation pressure of 10 ATM must be selected, and should be inflated to a minimum of 10 ATM.

## Controlled expansion

	Feature	Benefit
1	 ePTFE graft-lining	Sustained patency due to reduced in-stent stenosis
2	 Nitinol frame	Self-expanding, radial force and flexibility
3	 Unlined portal region	For portal perfusion and shunt perfusion
4	 Radiopaque gold markers	Easy visualization
5	 Controlled expansion sleeve	Size and set the diameter during implantation



### Time to rethink TIPS.

Learn more at [goremedical.com/rethinkTIPS](http://goremedical.com/rethinkTIPS)

Please contact your Gore technical sales associate if you have further questions.

 Consult Instructions for Use

Refer to *Instructions for Use* for a complete description of all warnings, precautions, and contraindications. <sup>®</sup> only

Products listed may not be available in all markets.

GORE, *Together, improving life*, VIATORR and designs are trademarks of W. L. Gore & Associates, Inc.  
© 2019 W. L. Gore & Associates, Inc. AUGUST 2019

W. L. Gore & Associates, Inc.  
Flagstaff, AZ 86004  
Asia Pacific +65 67332882 Europe 00800 6334 4673  
United States 800 437 8181 928 779 2771 [goremedical.com](http://goremedical.com)

